


Physcomitrella Measure Plugin

Find areas of high and low density in brightfield images of *Physcomitrella*

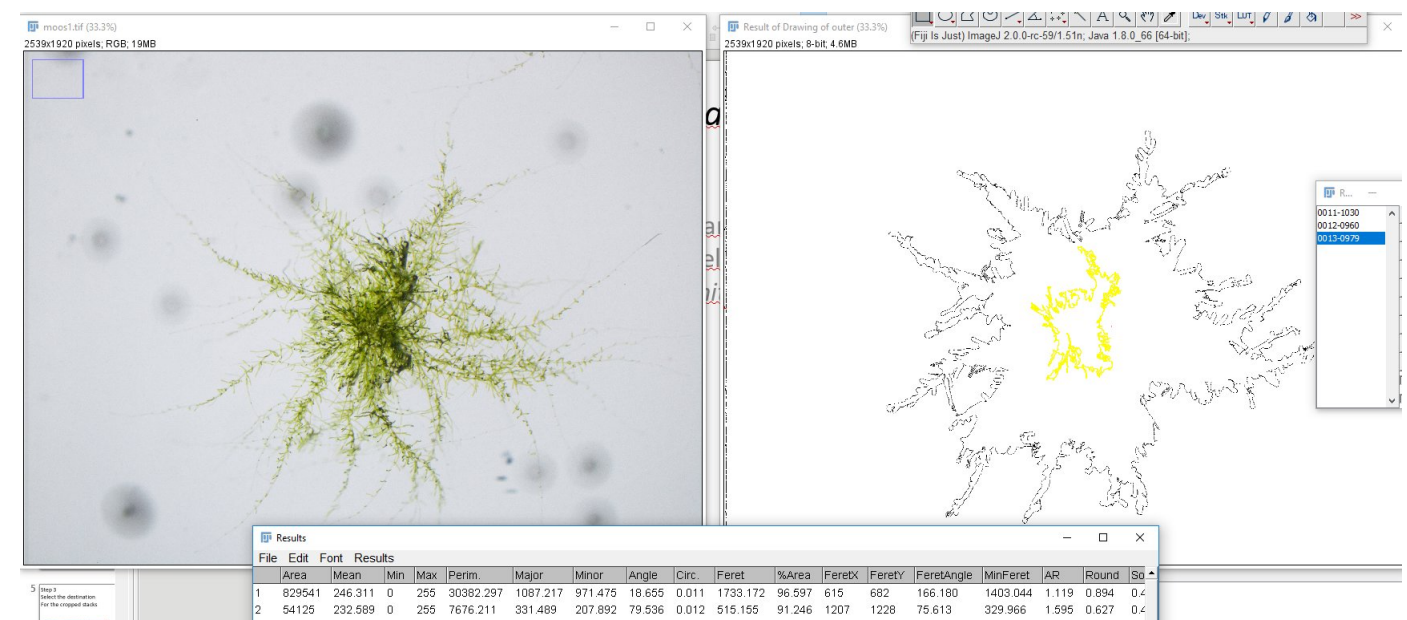

Simon Fraas

MLU Halle AG Quint

Simon.fraas@landw.uni-halle.de

Requirements

You need single brightfield images of *Physcomitrella* (or similar looking) material



20C_4



20C_5



20C_6



20C_7




20C_8




28C_1




28C_2

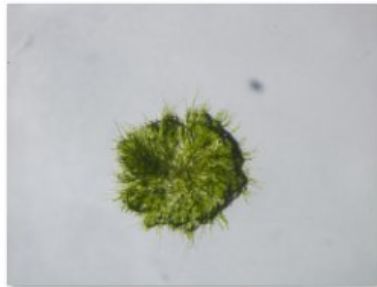

28C_3

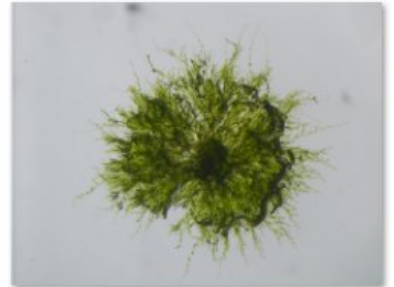


28C_4

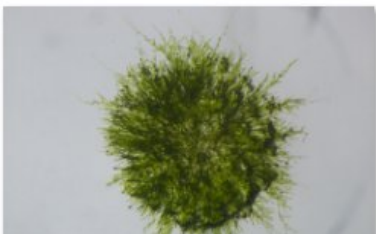


Step 1

Open picture

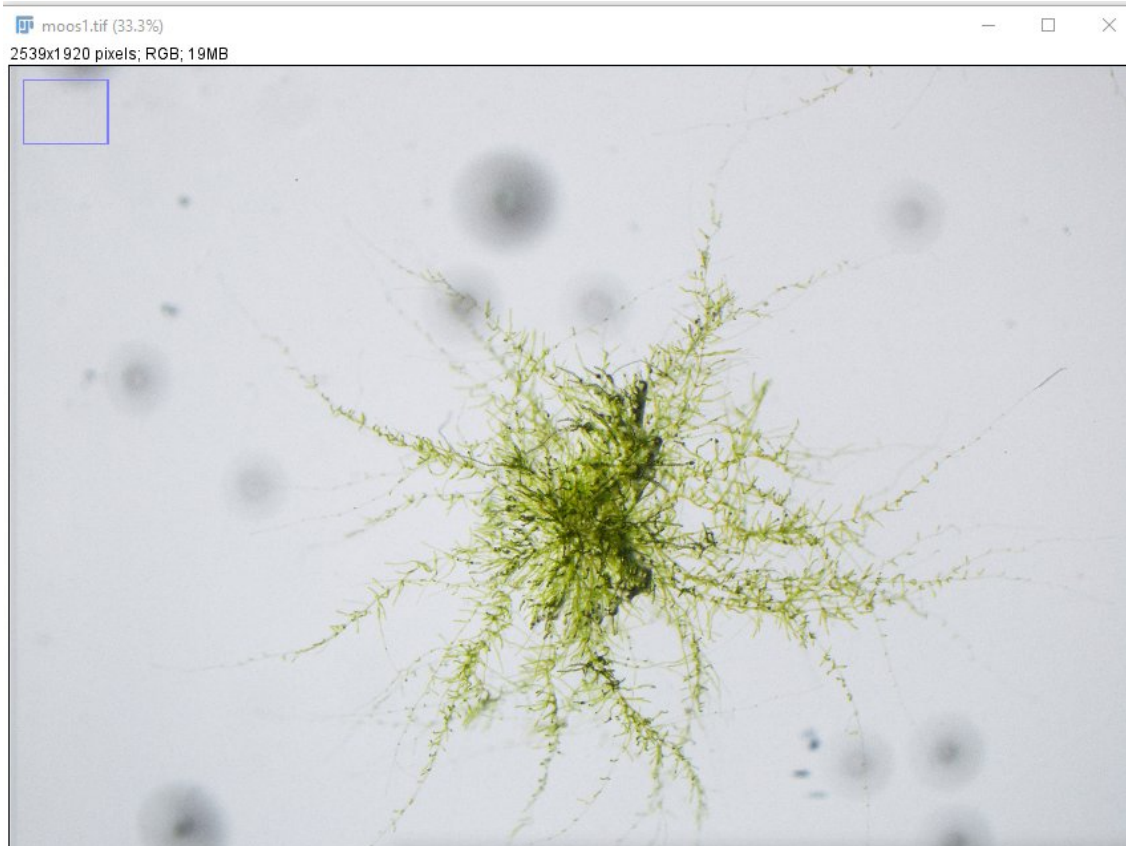

Step 2


Start the Plugin

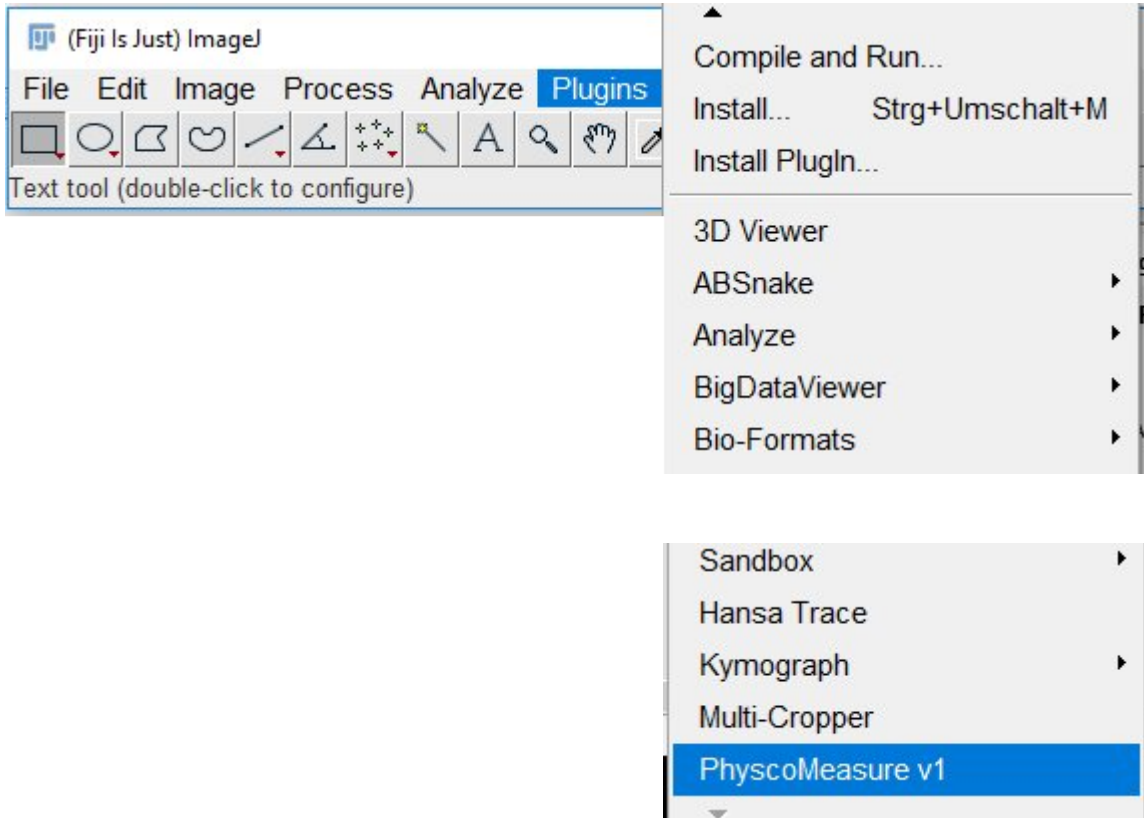

Step 3


Wait for the process to complete

In the end it should look like this

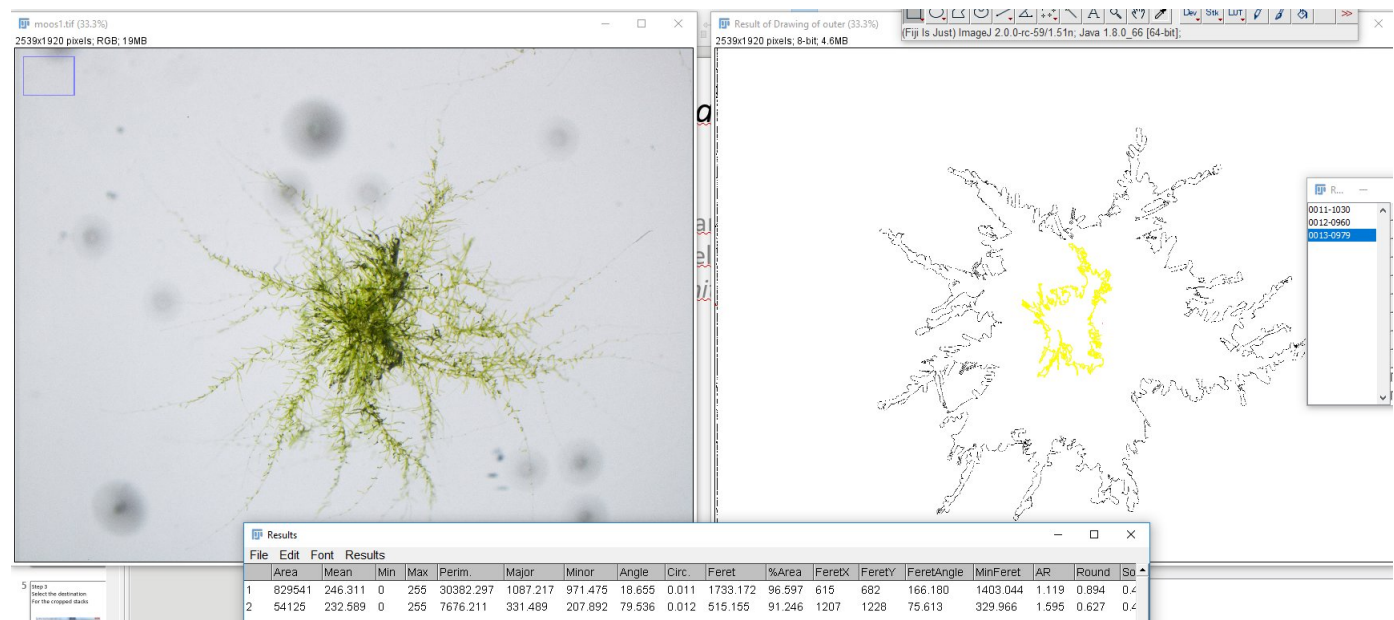

Select ROIs at the ROI manager window and hit M to measure.
Results can be saved as csv or copied directly

What actually happens ?

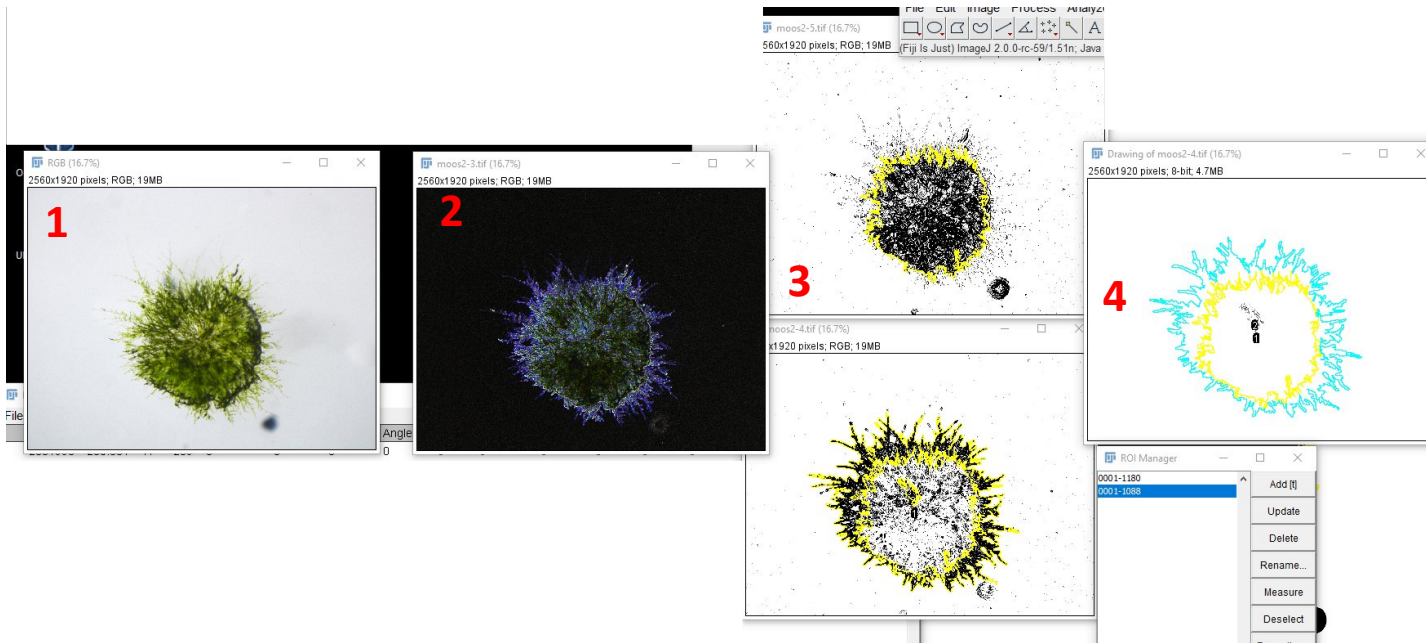

- 1 - A flat RGB will be created if you have a composite image.
- 2 - Edge detection via Sobel filter is applied
- 3 - Image is duplicated and color thresholds for dense and non dense area are run
- 4 - Thresholded areas are being extracted and measured via Particle filter
- 5 - User dialog view is being created

